

# TROPONIN I (HUMAN CARDIAC-SPECIFIC) ELISA

## Catalog Number: 1105



**DIAGNOSTIC AUTOMATION, INC.**

23961 Craftsman Road, Suite E/F,  
Calabasas, CA 91302

Tel: (818) 591-3030 Fax: (818) 591-8383

onestep@rapidtest.com

technicalsupport@rapidtest.com

www.rapidtest.com



See external label



2°C-8°C



Σ=96 tests



#1105

### Enzyme Immunoassay for the Quantitative Determination of Cardiac-Specific Troponin-I in Human Serum

#### FOR IN VITRO DIAGNOSTIC USE

Store at 2 to 8°C.

#### PROPRIETARY AND COMMON NAMES

Human Cardiac-Specific Troponin-I Enzyme Immunoassay  
(cTnI ELISA)

#### INTENDED USE

The cTnI ELISA is intended for the quantitative determination of cardiac troponin I in human serum. Measurement of troponin I values are useful in the evaluation of acute myocardial infarction (AMI).

#### SUMMARY AND EXPLANATION OF TEST

Troponin is the inhibitory or contractile regulating protein complex of striated muscle. It is located periodically along the thin filament of the muscle and consists of three distinct proteins: troponin I, troponin C, and troponin T.<sup>1-5</sup> Likewise, the troponin I subunit exists in three separate isoforms; two in fast-twitch and slow-twitch skeletal muscle fibers, and one in cardiac muscle.<sup>6-8</sup> The cardiac isoform (cTnI) is about 40% dissimilar, has a molecular weight of 22,500 daltons, and has 31 additional amino acid residues that are not present on the skeletal isoforms.<sup>3-4,8-12</sup> Antibodies made against this cardiac isoform are immunologically different from antibodies made against the other two skeletal isoforms,<sup>10,13</sup> and the unique isoform and tissue specificity of cardiac troponin I is the basis for its use as an aid in the diagnosis of acute myocardial infarction (AMI).<sup>2-4,8-9,13-17</sup>

Cardiac troponin I (cTnI) has been useful in the differential diagnosis of patients presenting to Emergency

Departments (ED) with chest pain.<sup>18-20</sup> Myocardial infarction is diagnosed when blood levels of sensitive and specific biomarkers, such as cardiac troponin, the MB isoenzyme of creatine kinase (CK-MB), and myoglobin, are increased in a clinical setting of acute ischemia.<sup>21-22,23</sup>

The most recently described and preferred biomarker for myocardial damage is cardiac troponin (I or T).<sup>23</sup> The cardiac troponins exhibit myocardial tissue specificity and high sensitivity. Likewise, cardiac TnI and CK-MB have similar release patterns (4-6 hours after the onset of pain), but the level of cTnI remains elevated for a much longer period of time (6-10 days), thus providing for a longer window of detection of cardiac injury.<sup>23-24</sup>

Normal levels of cTn I in the blood are very low. After the onset of an AMI, cTnI levels increase substantially and are measurable in serum within 4 to 6 hours, with peak concentrations reached in approximately 12 to 24 hours after infarction.<sup>24-28</sup> The

fact that cTnI remains elevated in serum for a much longer period of time, added to its enhanced diagnostic sensitivity and cardiac specificity, allows for the detection of AMI much earlier after the onset of ischemia (4 hours),<sup>1,25</sup> as well as the diagnosis of peri-operative infarction in situations where a high serum level of skeletal muscle proteins are expected.<sup>17</sup>

Additionally, recent data have identified a measurable relationship between cardiac troponin levels and long-term outcome after an episode of chest discomfort.<sup>24,29</sup> The studies suggest that the use of the cTnI demonstrates high predictive value in delineating the high risk group of unstable angina patients,<sup>30</sup> and that these tests may be particularly useful in evaluating patient condition prior to discharge from the ED.<sup>25,29,31</sup>

The cTnI Enzyme Immunoassay provides a rapid, sensitive, and reliable assay for the quantitative measurement of cardiac-specific troponin I. The antibodies developed for the test will determine a minimal concentration of 1.0 ng/ml, and there is no cross-reactivity with human cardiac troponin T or skeletal troponin T or I.

#### PRINCIPLE OF THE ASSAY

The cTnI ELISA test is based on the principle of a solid phase enzyme-linked immunosorbent assay. The assay system utilizes four unique monoclonal antibodies directed against distinct antigenic determinants on the molecule. Three mouse monoclonal anti-troponin I antibodies are used for solid phase immobilization (on the microtiter wells). The fourth antibody is in the antibody-enzyme (horseradish peroxidase) conjugate solution. The test sample is allowed to react simultaneously with the four antibodies, resulting in the troponin I molecules being sandwiched between the solid phase and enzyme-linked antibodies. After a 90-minute incubation at room temperature, the wells are washed with water to remove unbound-labeled antibodies. A solution of

tetramethylbenzidine (TMB) Reagent is added and incubated for 20 minutes, resulting in the development of a blue color. The color development is stopped with the addition of 1N hydrochloric acid (HCl) changing the color to yellow. The concentration of troponin I is directly proportional to the color intensity of the test sample. Absorbance is measured spectrophotometrically at 450 nm.

#### REAGENTS AND MATERIALS PROVIDED

1. Antibody-Coated Wells (1 plate, 96 wells)  
Microtiter wells coated with mouse monoclonal anti-TnI.
2. Reference Standard Set (1 set, 1.0 ml/vial)  
Contains 0, 2.0, 7.5, 30, and 75 ng/ml TnI, lyophilized.
3. cTnI Enzyme Conjugate Reagent (13 ml/vial)  
Contains mouse monoclonal anti-TnI conjugated to horseradish peroxidase in Tris Buffer-BSA solution with preservatives.
4. TMB Reagent (11 ml/bottle)  
Contains one-step TMB solution.
5. Stop Solution (11 ml/bottle)  
Contains diluted hydrochloric acid (1N HCl).

#### MATERIALS REQUIRED BUT NOT PROVIDED

1. Distilled or deionized water
2. Precision pipettes: 5  $\mu$ l, 10  $\mu$ l, 50  $\mu$ l, 100  $\mu$ l and 1.0 ml
3. Disposable pipette tips
4. Microtiter well reader capable of reading absorbance at 450 nm.
5. Vortex mixer, or equivalent
6. Absorbent paper
7. Graph paper
8. Cardiac Marker Tri-level Control; Cat. No. 685 (Bio-Rad Laboratories Diagnostics Group, Hercules, CA 94547)

#### WARNINGS AND PRECAUTIONS

1. CAUTION: This kit contains human material. The source material used for manufacture of this component tested negative for HBsAg, HIV 1/2 and HCV by FDA-approved methods. However, no method can completely assure absence of these agents. Therefore, all human blood products, including serum samples, should be considered potentially infectious. It is recommended that the reagents and patient samples be handled according to the OSHA Standard on Bloodborne Pathogens<sup>32</sup> or other appropriate national biohazard safety guidelines or regulations.<sup>33-34</sup>
2. Avoid contact with 1N HCl. It may cause skin irritation and burns. If contact occurs, wash with copious amounts of water and seek medical attention if irritation persists.
3. Do not use reagents after expiration date and do not mix or use components from kits with different lot numbers.
4. Replace caps on reagents immediately. Do not switch caps.
5. Do not pipette reagents by mouth.
6. For in vitro diagnostic use.

#### STORAGE CONDITIONS

1. Store the unopened kit at 2-8°C upon receipt and when it is not in use, until the expiration shown on the kit label. Refer to the package label for the expiration date.
2. Keep microtiter plate in a sealed bag with desiccant to minimize exposure to damp air.

#### REAGENT PREPARATION

1. All reagents should be allowed to reach room temperature (18-25°C) before use.
2. Reconstitute each lyophilized standard with 1.0 ml distilled water. Allow the reconstituted material to stand for at least 20 minutes and mix gently. ***The Reconstituted standards will be stable for up to 8 hours when stored sealed at 2-8 °C.*** Discard the reconstituted Standards after 8 hours. ***To assure maximum stability of the reconstituted Standards, they should be aliquoted and frozen (-20 °C or below) immediately after reconstitution has been achieved. Each aliquoted Standard should be frozen and thawed only once.***
3. Samples with expected Troponin I concentrations over 100 ng/ml may be quantitated by dilution with diluent available from vender.

#### INSTRUMENTATION

A microtiter well reader with a bandwidth of 10 nm or less and an optical density range of 0 to 3 OD or greater at 450 nm wavelength is acceptable for absorbance measurement.

#### SPECIMEN COLLECTION AND PREPARATION

1. The use of SERUM samples is required for this test.
2. Specimens should be collected using standard venipuncture techniques. Remove serum from the coagulated or packed cells *within 60 minutes* after collection.
3. Specimens which cannot be assayed within 24 hours of collection should be frozen at -20°C or lower, and will be stable for up to six months.
4. Avoid grossly hemolytic (bright red), lipemic (milky), or turbid samples (after centrifugation).
5. Specimens should not be repeatedly frozen and thawed prior to testing. DO NOT store in "frost free" freezers, which may cause occasional thawing. Specimens which have been frozen, and those which are turbid and/or contain particulate matter, must be centrifuged prior to use.

### PROCEDURAL NOTES

1. Pipetting Recommendations (single and multi-channel): Pipetting of all standards, samples, and controls should be completed within 3 minutes.
2. All standards, samples, and controls should be run in duplicate concurrently so that all conditions of testing are the same.
3. It is recommended that the wells be read within 15 minutes following addition of Stop Solution.

### ASSAY PROCEDURE

1. Secure the desired number of coated wells in the holder.
2. Dispense 100  $\mu$ l of standards, specimens, and controls into appropriate wells.
3. Dispense 100  $\mu$ l of Enzyme Conjugate Reagent into each well.
4. Thoroughly mix for 30 seconds. It is very important to mix completely.
5. Incubate at room temperature (18-25°C) for 90 minutes.
6. Remove the incubation mixture by flicking plate contents into a waste container.
7. Rinse and flick the microtiter wells 5 times with distilled or deionized water. (Please do not use tap water.)
8. Strike the wells sharply onto absorbent paper or paper towels to remove all residual water droplets.
9. Dispense 100  $\mu$ l of TMB Reagent into each well. Gently mix for 5 seconds.
10. Incubate at room temperature for 20 minutes.
11. Stop the reaction by adding 100  $\mu$ l of Stop Solution to each well.
12. Gently mix for 30 seconds. ***It is important to make sure that all the blue color changes to yellow color completely.***
13. Read absorbance at 450nm with a microtiter well reader ***within 15 minutes***

### QUALITY CONTROL

1. Good laboratory practice requires that quality control specimens (controls) be run with each calibration curve to verify assay performance. To ensure proper performance, control material should be assayed repeatedly to establish mean values and acceptable ranges.

### CALCULATION OF RESULTS

1. Calculate the mean absorbance value (OD<sub>450</sub>) for each set of reference standards, controls and samples.
2. Construct a standard curve by plotting the mean absorbance obtained for each reference standard against its concentration in ng/ml on graph paper, with absorbance on the vertical (y) axis and concentration on the horizontal (x) axis.
3. Using the mean absorbance value for each sample, determine the corresponding concentration of troponin I (ng/ml) from the standard curve. Depending on experience and/or the availability of computer capability, other methods of data reduction may be employed.


4. Patient samples with cTnI concentrations greater than 100 ng/ml should be diluted 10-fold with vender's Troponin I Sample Diluent. The final cTnI values should be multiplied by 10 to obtain cTnI results in ng/ml.

### EXAMPLE OF STANDARD CURVE

Results of a typical standard run with absorbency readings at 450nm shown on the Y axis against troponin I concentrations shown on the X axis. **NOTE:** This standard curve is for the purpose of illustration only, and should not be used to calculate unknowns. Each laboratory must provide its own data and standard curve in each experiment.

#### A. Standard Curve from Assay Procedure A:

cTnI (ng/ml)	Absorbance (450 nm)
0	0.048
2.0	0.110
7.5	0.307
30	1.357
75	2.853

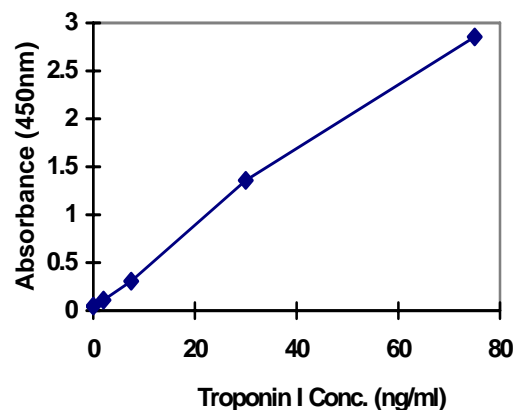

### LIMITATIONS OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the package insert instructions and with adherence to good laboratory practice.
2. Diagnostic results obtained from the cTnI ELISA should be used in conjunction with other diagnostic procedures and information available to the physician.; e.g., additional clinical testing, ECG, symptoms, and clinical observations.
3. Serum samples demonstrating gross lipemia, gross hemolysis, or turbidity should not be used with this test.
4. The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.
5. Patient samples may contain human anti-mouse antibodies (HAMA) which are capable of giving falsely elevated or depressed results with assays that utilize mouse monoclonal antibodies. The vender's cTnI

ELISA assay has been designed to minimize interference from HAMA-containing specimens; nevertheless complete elimination of this interference from all patient specimens cannot be guaranteed.

6. Test results that are inconsistent with the clinical picture and patient history should be interpreted with caution.

### EXPECTED VALUES

An evaluation of the clinical data was conducted to determine the normal expected value, as well as the clinical sensitivity and clinical specificity of the vebder's cTnI assay (see below). Two-hundred and twenty-five (225) apparently healthy adults were assayed using the test to establish the normal expected value, which was determined to be  $\leq 0.5$  ng/ml cTnI. All values from the normal population tested were below the sensitivity level of the assay (1.0 ng/ml).

It is recommended that each laboratory establish its own normal range based on the patient population, geography, dietary and environmental factors; likewise current practice and clinical criteria for AMI diagnosis must be considered. However, based on published literature, the diagnostic cut-off for AMI patients is determined to be 1.5 ng/ml.<sup>35</sup>

Any conditions resulting in myocardial cell damage can potentially increase cardiac troponin-I levels above the expected value. These conditions have been documented clinically to include unstable angina, myocarditis, congestive heart failure, and cardiac surgery or invasive testing.<sup>3,29</sup>

**NOTE:** Serial sampling may be required to detect elevated levels.

### CLINICAL PERFORMANCE

A clinical investigation was conducted to determine the accuracy, as well as the diagnostic sensitivity and specificity of the cTnI ELISA as compared to another commercially available kit. The data is presented below.

#### 1. Clinical Correlation

A statistical study using 204 clinical patient serum samples, ranging in cTnI concentration from 0.7 ng/ml to 595 ng/ml as analyzed using the cTnI ELISA (0.5 ng/ml to 484 ng/ml; Abbott TnI MEIA), demonstrated equivalent correlation with a commercially available kit as shown below.

Comparison between the cTnI ELISA and the Abbott AxSym<sup>®</sup> TnI test provided the following data:

Correlation coefficient = 0.9537  
 Slope = 0.9063  
 Intercept = -3.9875  
 Mean = 45.57 ng/ml  
 Abbott Mean = 50.37 ng/ml

When samples which were above the upper limit of the Abbott assay were removed (i.e., > 50 ng/ml), the following statistics were observed. Note this was done

in order to demonstrate the concordance between assays in an undiluted sample population:

Correlation coefficient = 0.8672  
 Slope = 1.0416  
 Intercept = 0.7816  
 Mean = 9.96 ng/ml  
 Abbott Mean = 12.72 ng/ml

#### 2. Clinical Sensitivity and Specificity

Of the overall 149 total patients (249 samples) evaluated in the study, there were 93 patients who were confirmed to have experienced an AMI. Based on a clinical cutoff of 1.5 ng/ml, the diagnostic sensitivity and specificity of the cTnI assay were evaluated. Clinical specificity was reported as 87.5% (95%CI: 80.3% - 94.7%), while sensitivity was calculated at 100%.

The results of these investigations demonstrate that the cTnI ELISA has comparable diagnostic accuracy to that of another currently marketed device.

### PERFORMANCE CHARACTERISTICS

#### 1. Sensitivity

The minimum detectable concentration of the cTnI ELISA assay as measured by 2SD from the mean of a zero standard is estimated to be 1.0 ng/ml. Additionally, the functional sensitivity was determined to be 0.75 ng/ml ((as determined with inter-assay %C.V.  $\leq 10\%$ ). **Lower limit of cTnI ELISA  $\cong 0.48$  ng/ml cTnI; upper limit = 1.0 ng/ml cTnI.**

#### 2. Hook Effect

No hook effect was observed in this assay at cardiac troponin-I concentrations up to 10,000 ng/ml.

#### 3. Precision

##### a. Intra-Assay Precision

Within-run precision was determined by replicate determinations of seven different serum samples in one assay. Within-assay variability is shown below:

Serum Sample	1	2	3	4	5
# Repts.	20	20	20	20	20
Mean cTnI (ng/ml)	1.69	5.93	24.3	44.9	89.8
S.D.	0.05	0.22	1.35	1.78	2.52
C.V. (%)	2.8	3.7	5.6	4.0	2.8

##### b. Inter-Assay Precision

Between-run precision was determined by replicate measurements of six different serum samples over a series of individually calibrated assays. Between-assay variability is shown below:

Serum Sample	1	2	3	4	5
# Replicates	20	26	26	26	26
Mean cTnI (ng/ml)	1.68	5.88	24.56	48.91	85.81
S.D.	0.11	0.28	1.14	2.23	3.76
C.V. (%)	6.7	4.8	4.7	4.6	4.4

#### 4. Specificity

The following were tested for cross-reactivity at concentrations up to the levels indicated below. No cross-reactivity was observed for any of the components.

MATERIAL TESTED	TEST CONCENTRATION
Rabbit skeletal muscle troponin C	2,500 ng/ml
Human cardiac troponin T	2,500 ng/ml
Human skeletal muscle troponin T	2,500 ng/ml
Human skeletal muscle troponin I	2,500 ng/ml
Hemoglobin	1.2 g/dl
Bilirubin	20 mg/dl
Cholesterol	500 mg/dl
Triglyceride	1,000 mg/dl
Total Protein	10 g/dl

#### 5. Recovery and Linearity Studies

##### a. Recovery

Various patient serum samples of known human cTnI levels were combined and assayed in duplicate. The mean recovery was 93.3%.

PAIR NO.	EXPECTED [cTnI] (ng/ml)	OBSERVED [cTnI] (ng/ml)	% RECOVERY
1	4.25	3.95	92.9%
2	8.97	8.50	94.8%
3	11.43	10.49	91.8%
4	14.97	14.00	93.5%
5	32.34	29.62	91.6%
6	32.77	30.49	93.0%
7	81.00	77.21	95.3%

#### PERFORMANCE CHARACTERISTICS

##### b. Linearity

Four patient samples were serially diluted to determine linearity. The mean recovery was 101.7%.

#	Dilution	Expected Conc. (ng/ml)	Observed Conc. (ng/ml)	% Expected
1.	Undiluted	-----	-----	-----
	1:2	74.9	74.9	100%
	1:4	37.5	37.4	99.7%
	1:8	18.7	19.3	103.2%
	1:16	9.4	9.9	105.3%
	1:32	4.7	5.0	106.4%
	1:64	2.4	2.6	108.3%
	1:128	1.2	1.3	108.3%
Mean = 105.3%				
2.	Undiluted	-----	-----	-----
	1:2	68.4	68.4	100.0%
	1:4	34.2	34.2	100.0%
	1:8	17.1	17.6	102.9%
	1:16	8.6	8.5	98.8%
	1:32	4.3	4.4	102.3%
	1:64	2.2	2.4	109.1%
	1:128	1.1	1.2	109.1%
Mean = 103.2%				
3.	Undiluted	-----	-----	-----
	1:2	-----	-----	-----

	1:4	62.5	64.0	102.4%
	1:8	31.3	31.3	100.0%
	1:16	15.6	14.5	92.9%
	1:32	7.8	7.2	92.3%
	1:64	3.9	3.7	94.9%
	1:128	1.9	2.0	105.3%
Mean = 98.0%				
4.	Undiluted	-----	-----	-----
	1:2	-----	-----	-----
	1:4	86.4	88.0	101.9%
	1:8	43.2	43.1	99.8%
	1:16	21.6	21.7	100.5%
	1:32	10.8	10.2	94.4%
	1:64	5.4	5.4	100.0%
	1:128	2.7	2.8	103.7%
Mean = 100.1%				

#### STANDARDIZATION

Human Troponin I-T-C Complex was obtained from a qualified vendor, and cTnI concentration was determined. The material was further diluted with the cTnI Sample Diluent and served as "Standard Stock Solution" for preparing cTnI reference Standard Sets. The target value of the "Standard Stock Solution" was confirmed by the Abbott AxSym Troponin I immunoassay.

#### REFERENCES

- Etievent, J., Chocron, S., Toubin, G., et. al.: The use of cardiac troponin I as a marker of peri-operative myocardial ischemia. *Ann. Thorac. Surg.*, 59:1192-94, 1995.
- Apple, F.: Acute Myocardial Infarction and Coronary Reperfusion: Serum Cardiac Markers for the 1990's. *Am. J. Clin. Path.*, 97: 217-26, 1992.
- Adams, J., Bodor, G., Davila-Romain, V., et. al.: Cardiac troponin I: A marker for cardiac injury. *Circulation*, 88:101-6, 1993.
- Corin S., Juhasz, O., Zhul, et. al.: Structure and expression of the human slow twitch skeletal muscle troponin I gene. *J. Bio. Chem.*, 269:10651-7, 1994.
- Perry, S.V.: The regulation of contractile activity in Muscle. *Biochem. Soc. Trans.* 7:593, 1979.
- William, J.M., and Grand, R.J.A.: Comparison of amino acid sequence of troponin I from different striated muscles. *Nature*, 271:31, 1978.
- Mehegan, J.P., and Tobacman, L.S.: Cooperative interaction between troponin molecules bound to the cardiac thin filament. *J. Biol. Chem.*, 266:966, 1991.
- Mair, J., Wagner, I., Puschendorf, B., et. al.: Cardiac troponin I to diagnose myocardial injury. *The Lancet*, 341:838-9, 1993.
- Mair, J., Laruc, C., Mair, P., et al.: Use of Cardiac Troponin I to Diagnose Perioperative Myocardial Infarction in Coronary Artery Bypass Grafting. *Clin. Chem.*, 40: 2066-70, 1994.
- Vallins, W.J., et al.: Molecular cloning of human cardiac troponin I using polymerase chain reaction. *FEBS Lett.*, 270:57, 1990.
- Leszyk, J., Dumaswala, R., Potter, J.D., et. al.: Amino acid sequence of bovine cardiac troponin I. *Biochemistry*, 27: 2821-7, 1988.
- Hartner, K.T., and Pette, D.: Fast and slow isoforms of troponin I and troponin C. Distribution in normal rabbit muscles and effects of chronic stimulation. *Eur. J. Biochem.* 188:261, 1990.
- Ebashi, S.: Ca<sup>2+</sup> and the contractile proteins. *J. Mol. Cell. Cardiol.* 16:129, 1984.
- Bodor, S.G., et. al.: Development of monoclonal antibody for an assay of cardiac troponin and T. Preliminary results in suspected cases of myocardial infarction. *Clin. Chem.* 38:2203, 1992.

15. Adams, J.E., et. al.: Biochemical markers of myocardial injury. Is MB creatine kinase the choice for the 1990's. *Circulation*, 88:750, 1993.
16. Hamm, C.W.: Cardiac-specific Troponins in Acute Coronary Syndromes, in Heart Disease, a textbook of Cardiovascular Medicine, ed. Braunwald, E., W.B. Saunders Co. Update 3, 1977.
17. Adams, J., Sicard, G., Allen, B., et. al.: Diagnosis of perioperative myocardial infarction with measurement of cardiac troponin I. *N. Eng. J. Med.*, 330: 670-4, 1994.
18. Meyer, T., Binder, L., Graeber, T., et. al.: Superiority of combined CK-MB and troponin I measurements for the early risk stratification of unselected patients presenting with acute chest pain. *Cardiology*, 90:286-94, 1998.
19. Polanczyk, C.A., Lee, T.H., Cook, E.F., et. al.: Cardiac troponin I as a predictor of major cardiac events in emergency department patients with acute chest pain. *J. Am. Coll. Cardiol.*, 32:8-14, 1998.
20. D'Costa, M., Fleming, E., Patterson, M.C.: Cardiac troponin I for the diagnosis of acute myocardial infarction in the emergency department, *Am. J. Clin. Pathol.*, 108:550-5, 1997.
21. Rice, M.S., MacDonald, D.C.: Appropriate roles of cardiac troponins in evaluating patients with chest pain. *J. Am. Board Fam. Pract.*, 12:214-8, 1999.
22. Hamm, C.W., Goldmann, B.U., Heeschen, C., et. al.: Emergency room triage of patients with acute chest pain by means of rapid testing for cardiac troponin T or troponin I. *N. Engl. J. Med.*, 337:1648-53, 1997.
23. Joint European Society of Cardiology/American College of Cardiology: *J. Am. Coll. Cardio.*, 36(3), "Myocardial Infarction Redefined, 2000.
24. Cummins, B., et. al.: Cardiac-specific troponin-I radioimmunoassay in the diagnosis of acute myocardial infarction. *Am. Heart J.*, 113:1333, 1987.
25. BLOOD TESTS FOR RAPID DETECTION OF HEART ATTACK © 2000 American Heart Association, Inc. Guideline, September, 2000.
26. Ebell, M.H., White, L.L., and Weismantel, D.: A systematic review of troponin T and I values as a prognostic tool for patients with chest pain. *J. Fam. Pract.*, 49:746-53, 2000.
27. Mair, J., Morandell, D., Gonser, N., et. al.: Equivalent early sensitivities of myoglobin, creatine kinase-MB mass, creatine kinase, isoform rallo, and cardiac troponins I and T for acute myocardial infarction. *Clin. Chem.*, 41: 1266-72, 1995.
28. Bertinchant, J.P., Laruc, C., Pemel, I., et. al.: Release kinetics of serum cardiac troponin I in ischemic myocardial injury. *Clin. Bio. Chem.*, 29:587-94, 1996.
29. Bodor, G.S.: Cardiac troponin I: a highly biochemical marker for myocardial infarction. *J. Clin. Immunoassay*, 17: 40-4, 1994.
30. Antman, E.M., Tanasijevic, M.J., Thompson, B., et. al.: Cardiac-specific troponin I levels to predict the risk of mortality in patients with acute coronary syndromes. *N. Engl. J. Med.*, 335:1342-9, 1996.
31. Ellestad, M.H., The diagnostic power of four chemical markers on admission to the chest pain center. In *Differential Diagnosis and Management of Patients with Chest Pain: A Multiple Biochemical Marker Approach*. A Symposium Sponsored by Baylor College of Medicine, held during the 68th Scientific Session of the American Heart Association, November 11, 1995, Anaheim, California.
32. Hanfner, S., et. al.: Cardiac troponins in serum in chronic renal failure. *Clin. Chem.*, 40:1790 1994.
33. U.S. Department of Labor, Occupational Safety and Health Administration, 29 CFR Part 1910.1030. Occupational Exposure to Bloodborne Pathogens; Final Rule. *Federal Register*; 56(235):64175, 1991.
34. USA Center for Disease Control/National Institute of Health Manual, "Biosafety in Microbiological and Biomedical Laboratories", (1984).
35. National Committee for Clinical Laboratory Standards. Protection of Laboratory Workers from Instrument Biohazards and Infectious Disease Transmitted by Blood, Body Fluids, and Tissue: Approved Guideline. NCCLS Document M29-A, 1997.
36. "Clinical Guide to Laboratory Tests". Edited by N.W. Tietz, 3rd Edition. W.B. Saunders Company, Philadelphia, PA. 19106 (1995).



**DIAGNOSTIC AUTOMATION, INC.**  
 23961 Craftsman Road, Suite E/F, Calabasas, CA 91302  
 Tel: (818) 591-3030 Fax: (818) 591-8383  
**ISO 13485-2003**

Revision Date: 4/6/06